




## SEPAX APPLICATION NOTE

# Nucleotide Analysis on Proteomix SAX Column

**Proteomix SAX, NP5, 4.6 x 100 mm**  
Part Number: 403NP5-4610

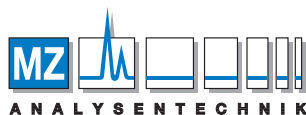

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: [info@mz-at.de](mailto:info@mz-at.de), [www.mz-at.de](http://www.mz-at.de)



# Nucleotide Analysis on Proteomix SAXNP5 Column

Column: Proteomix SAX, NP5, 4.6 x 100 mm; Instrument: Agilent 1260; Mobile Phase: A: 20 mM Tris, pH 8, B: A+ 1 M NaCl; Gradient: 0-20% B in 15 min; Flow Rate: 0.5 mL/min; Pressure: 72 bar; Injection Volume: 2  $\mu$ L; Detection: 260 nm; Column Temperature: 23  $^{\circ}$ C; Pressure: 72 bar

