


Insecticides (Neonicotinoids) Analysis (DE-413)


Neonicotinoids, one category of insecticides, were analyzed using RSpak DE-413 .

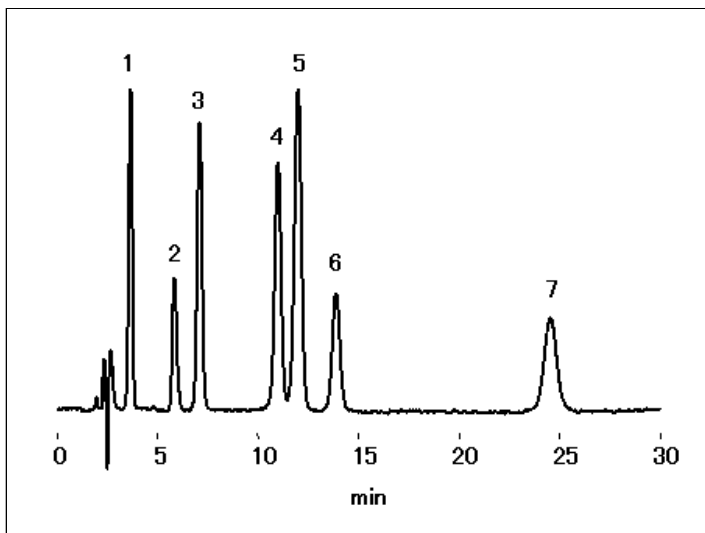

Sample : 1.25µg/mL each, 10µL

1. Dinotefuran
2. Nitenpyram
3. Thiamethoxam
4. Clothianidin
5. Acetamiprid
6. Imidacloprid
7. Thiacloprid

Column : Shodex RSpak DE-413 (4.6mmI.D. x 150mm)
Eluent : 10mM Formic acid aq./CH₃CN=80/20
Flow rate : 0.8mL/min
Detector : Shodex RI
Column temp. : 40°C

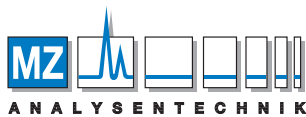

AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0, Fax +49 6131 880 96-20
e-mail: info@mz-at.de, www.mz-at.de