


HPLC Separation of perfluoroalkyl acids PFBA, PFAO, PEAS

Sample 21.7.2020 Owned by method@sielc.com

~ 2

Sample injection 21.7.2020 21:36

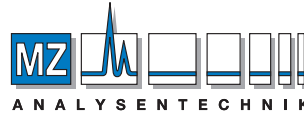

MZ-Analysentechnik GmbH
Barcelona-Allee 17 • D-55129 Mainz
Tel +49 6131 880 96-0
Fax +49 6131 880 96-20
e-mail: info@mz-at.de
www.mz-at.de

AUTHORIZED DISTRIBUTOR

Description

Imported from PFAO.CSV

Method 17 min, 1.0000 µl, 1.0000 mL/min, MeCN/H2O - 40/60%


Run time **17 min** Injection **1.0 µl** in 0.1 sec
Mobile phase ... MeCN/H2O - 40/60% CAD (Analog signal)
Buffer AmFm pH3.0 – 50mM
Flow rate **1.00 mL/min**

Column Newcrom A NA-46.150.0510

Serial NAQ8QWXA Length × ID 150 × 4.6 mm
Type NA-46.150.0510 Particle, pore size 5 µm, 100 Å
Lot B04-051 Max pressure 500 psi

Data

PFAO ...

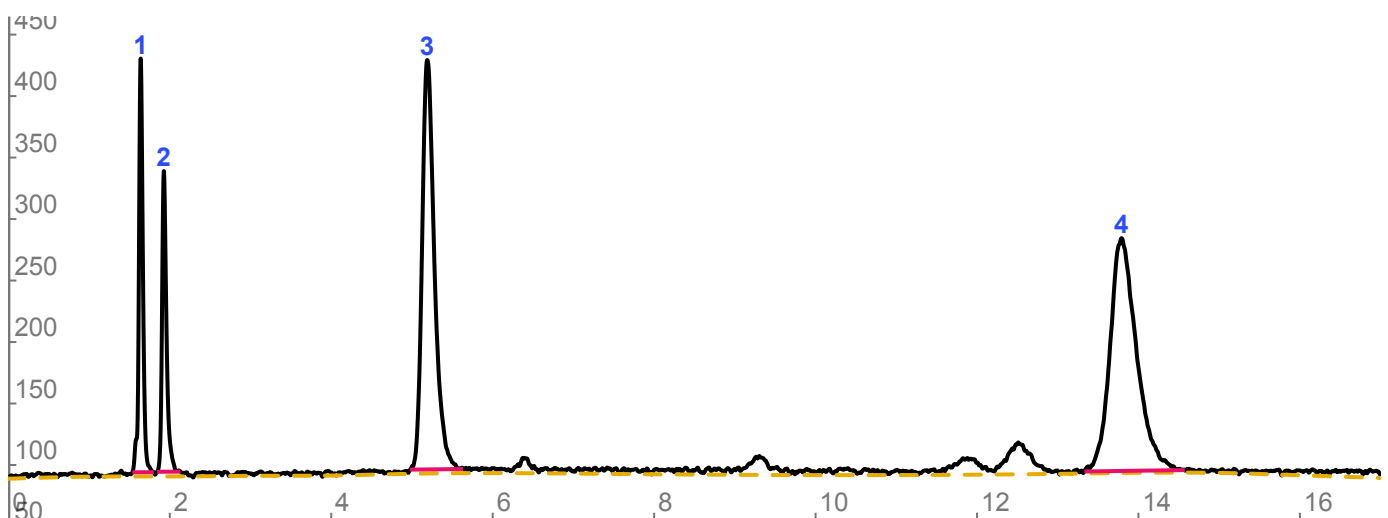

#	Time	Area	Compound	Amount
1	1.64	19.442	→ K	
2	1.93	15.741	→ Perfluorobutyric acid	