


DATA CODE:

T01130

Application Title:  
Separation of DNA fragments

Compound:  
1kb DNA Ladder

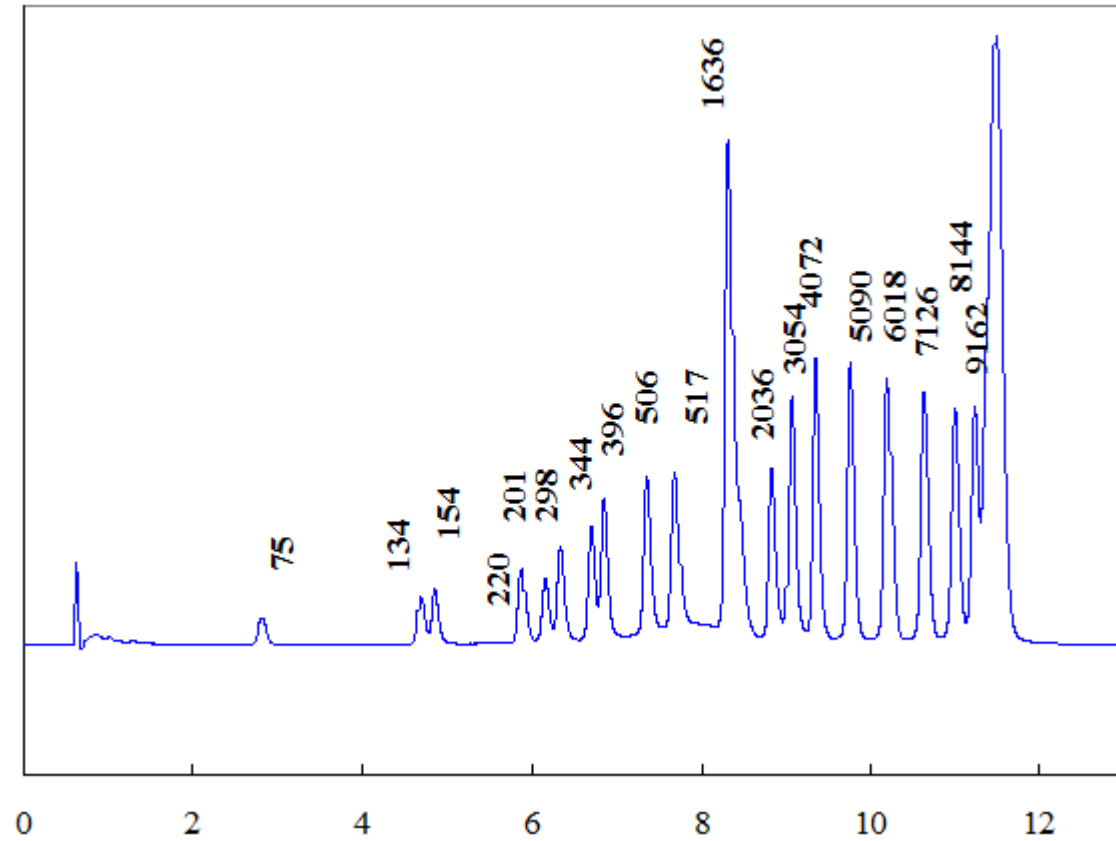

### Chromatographic conditions

Column:	TSKgelDNA-NPR 4.6mm×7.5cm, P/N 18249, 4.6mm×7.5cm
Mobile Phase:	A:20mmol/Ltris-HCl (pH9.0) ; B:20mmol/Ltris-HCl (pH9.0) +1mol/L NaCl; gradient:0min (B50%)→20min (B75%)
Flow Rate:	0.75mL/min
Detection:	UV/VIS 260nm
Temperature:	25°C
Injection Volume:	10µL
Title:	Separation of DNA fragments