

Data code:


G00722

Application Title:

Separation of low molecular mass heparin by size-exclusion chromatography.

Compound:

low molecular mass heparin


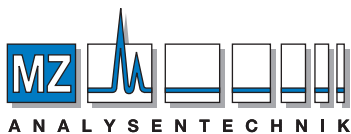


AUTHORIZED DISTRIBUTOR

MZ-Analysentechnik GmbH, Barcelona-Allee 17• D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de

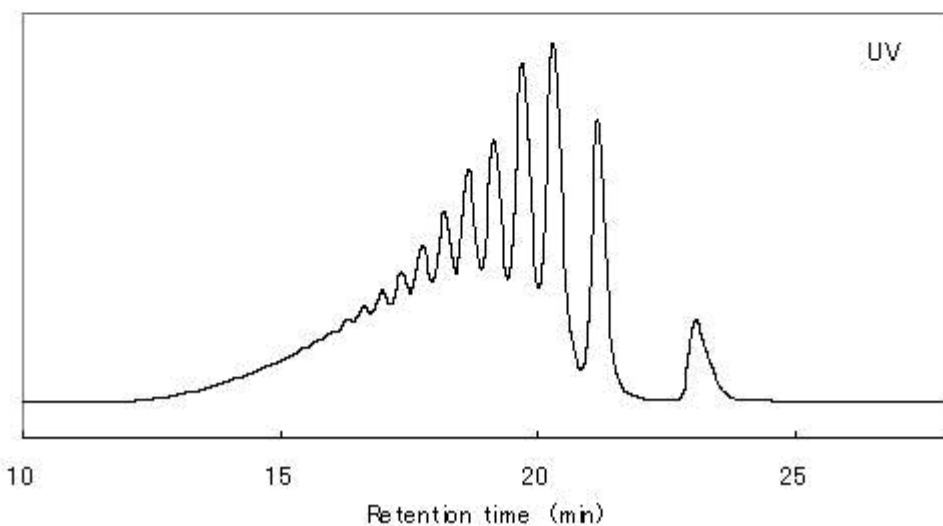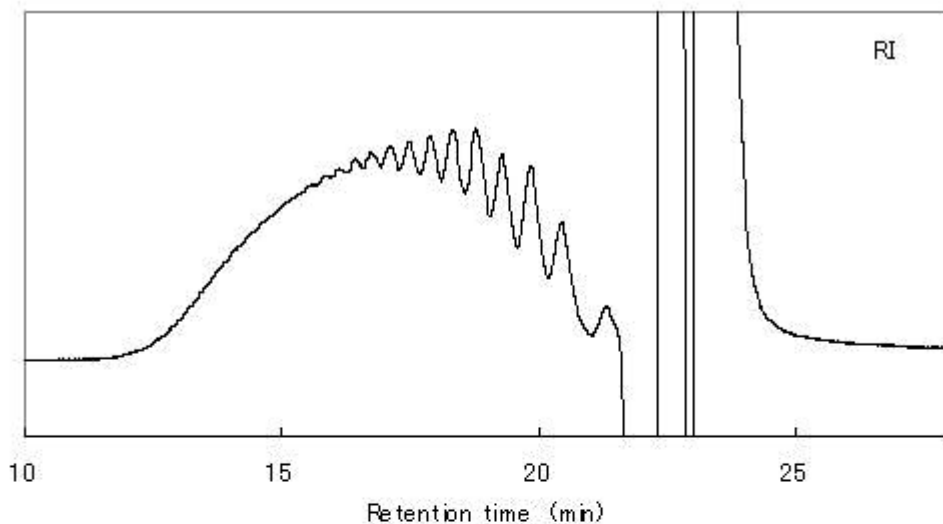

Chromatographic conditions

Column:	TSKgel G2000SW _{XL} , P/N 08540, 5 μ m, 7.8 mm ID \times 30 cm \times 1
Mobile Phase:	0.2 mol/L sodium sulfate (pH 5.0)
Flow Rate:	0.5 mL/min
Detection:	RI, UV/VIS @ 234 nm
Temperature:	25 $^{\circ}$ C
Injection Volume:	50 μ L
Sample(s):	Standard
Sample Load:	1 g/L