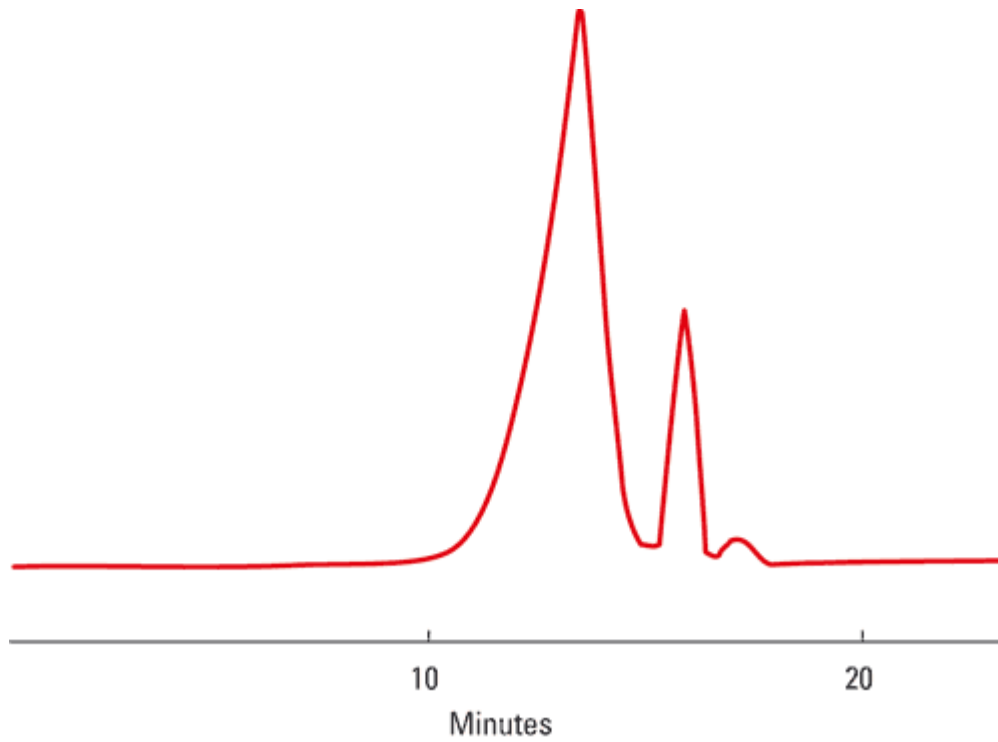


Application Title:Chromatogram of DNA-Na by-size exclusion chromatography on TSKgel GMPW_{XL}**Compound:**

DNA-Na

**Chromatographic conditions**

Column:	TSKgel GMPW _{XL} , P/N 08025, 13 µm, 7.8 mm ID × 30 cm × 2
Mobile Phase:	0.2 mol/L phosphate buffer (pH 6.9)
Flow Rate:	1 mL/min
Detection:	RI
Sample(s):	Polymer
Reference:	TOSOH Separation Report
Title:	SEC analysis of water soluble polymer for TSKgel PW type column (2) -PWXL series-
Number:	38
Page start:	12
Page end:	12