

Data code:


T01015

Application Title:

Separation of acrylamide by RPC

Compound:

acrylamide; acrylic acid

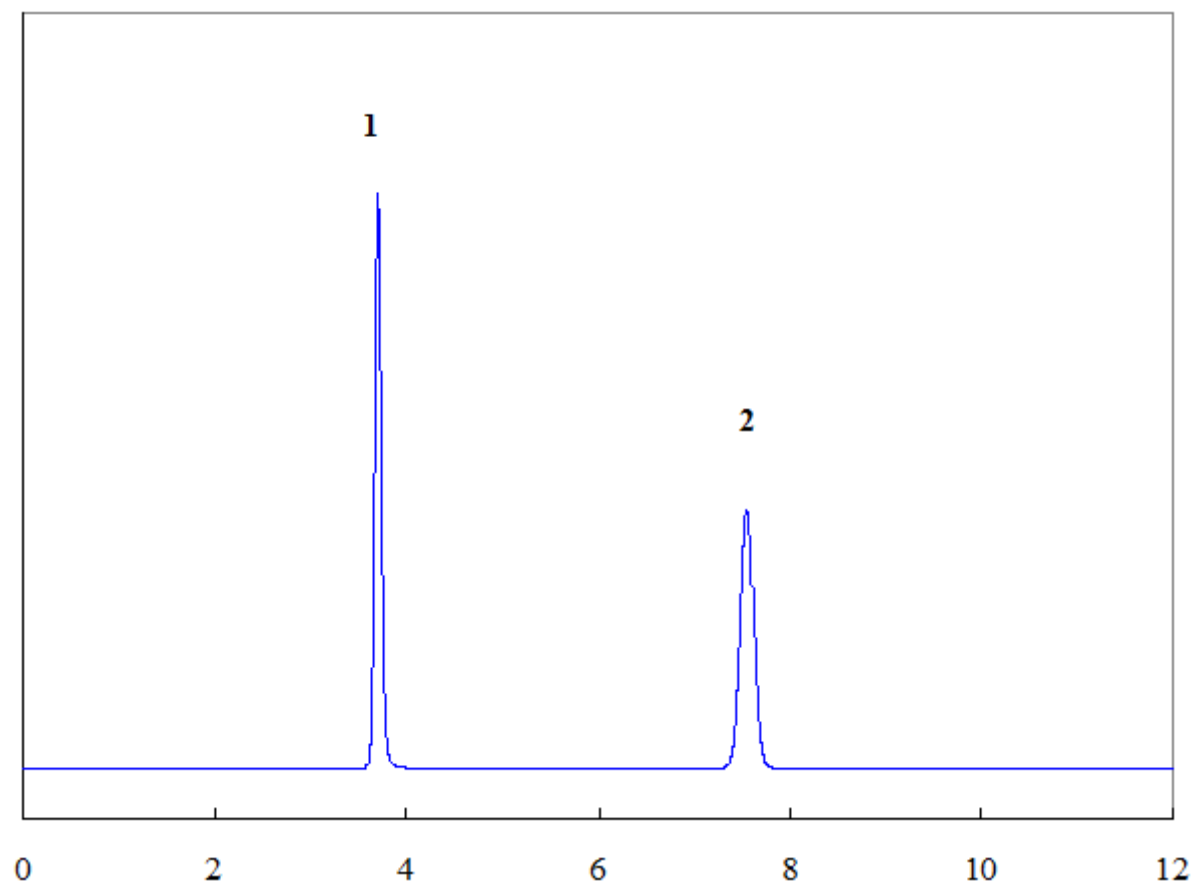

Chromatographic conditions

Column: TSKgelODS-100V 5 μ m 4.6mm \times 15cm, P/N 21455, 4.6mm \times 15cm

Mobile Phase: 0.1%H₃PO₄

Flow Rate: 1.0mL/min

Detection: UV/VIS 210nm

Temperature: 40 $^{\circ}$ C

Injection Volume: 5 μ L

Title: Separation of acrylamide by RPC