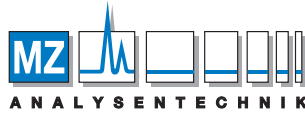


Application Title:

Separation of DNA from clupeoid seminal fluid

Compound:


DNA from clupeoid seminal fluid

**AUTHORIZED DISTRIBUTOR**

MZ-Analysentechnik GmbH, Barcelona-Allee 17 • D-55129 Mainz

Tel +49 6131 880 96-0, Fax +49 6131 880 96-20

e-mail: info@mz-at.de, www.mz-at.de

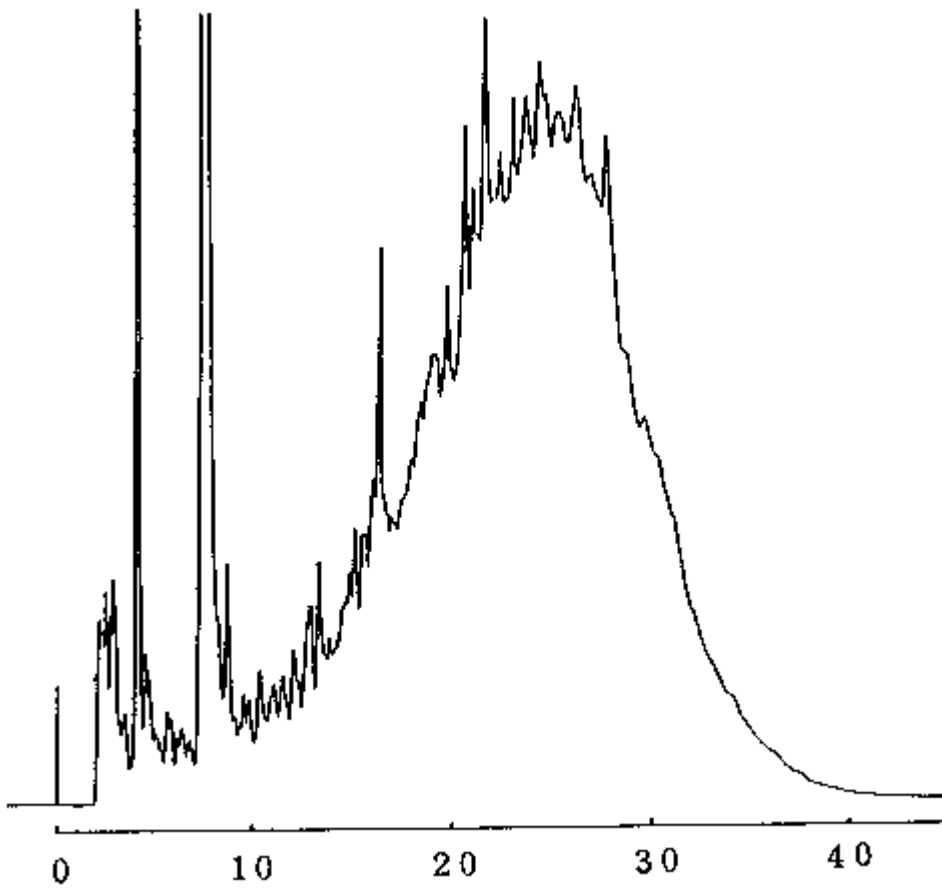**Chromatographic conditions**

Column:	TSKgel OligoDNA RP, P/N 13352, 5 µm, 4.6 mm ID × 15 cm × 1
Mobile Phase:	A: Acetonitrile / 0.1 mol/L ammonium acetate buffer (pH 7.0) = 2.5 / 97.5 B: Acetonitrile / 0.1 mol/L ammonium acetate buffer (pH 7.0) = 25 / 75 Linear gradient from A to B for 90min
Flow Rate:	1 mL/min
Detection:	UV/VIS @ 260 nm
Sample(s):	Standard
Reference:	TOSOH Separation Report
Title:	Separation of nucleic acid on TSKgel Oligo-DNA RP
Number:	51
Page start:	8

